

NEW! Bathmaster Sonaris™

An economical reclining bath lift that provides comfort, safety and independence.
Includes improved warranty and new accessories.

Being able to bathe independently while feeling safe and comfortable is now more affordable than ever. The new Sonaris is the latest in a line of popular Bathmaster bath lifts that delivers a superior bathing experience with features such as:

- Battery-powered motor so no installation is required
- Seat and back simply click together for quick assembly, or easily detach for portability
- Raising and lowering action is activated effortlessly by pressing the up or down pads on the hand controller
- Users can recline the backrest for greater submersion of their upper body
- Easy-clean, all-plastic design is lightweight with the heaviest piece weighing only $\frac{3}{4}$ lbs.
- 5-year warranty on frame. 2-year warranty on hand controller, recharger and cover.



Shown with the standard cover, 5670-10



Infinitely adjustable to 40°



Up to 18 lifts per charge

Here's a few more reasons why users enjoy the Sonaris:

- They can select any angle of recline up to 40° to achieve the most comfortable and safe position
- Recline does not occur automatically as it does on other lifts, but is entirely controlled by the user after the lift has completed its descent
- Speed is decreased during the reclining of the backrest so the user feels more comfortable
- A built-in safety mechanism will not allow users to lower the bath lift if the battery is too low to raise them
- An audible beep and warning light identifies when the battery requires charging and when to stop using the bath lift
- Superior battery provides an amazing 18 lifts per charge (based on 220-lb. load)
- It only requires 3 hours to fully charge, so it'll be ready when the user is
- Holding points on each component simplify carrying and transportation

5670-07 Sonaris Bathlift. Includes a hand controller with battery and recharger. Optional covers are sold separately. **Hy**

5670-14 Sonaris Replacement Hand Controller. Incorporates a superior battery and audible/visual alarms to signify operation and battery charge levels. Also has a handy sleep function to preserve the battery life. **Hy**

5670-15 Sonaris Replacement Recharger. Plug recharger into a standard outlet and the hand controller will be fully charged in about three hours. **Hy**

5670-10 Sonaris Standard Cover. Attractive ¼" closed-cell foam pad snaps into place on both the seat and the back providing additional comfort to the bather. **Hy**



NEW Accessories. **Hy**

5642-72 Swivel seat (sold without seat cover)

5642-73 Pair of replacement transfer flaps


5642-74 Set of 4 replacement feet

5642-70 Headrest (sold without headrest cover)

5642-71 Safety belt

5652-22 Headrest Blue Cover

5652-23 Swivel Seat Blue Cover

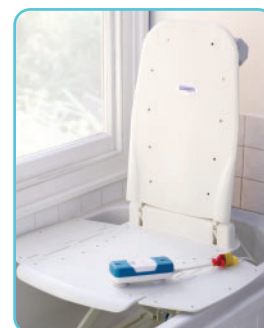

Bathmaster Sonaris Specifications

Seat Size:	20" x 27" (including side flaps)
	20" x 15" (excluding flaps)
Backrest Size:	24" x 14"
Height Range:	3"-18"
Seat Weight:	15¼ lbs.
Backrest Weight:	10¼ lbs.
Maximum User Weight:	309 lbs.
Backrest Reclining Speed:	3.93°/sec.
Average # of lift cycles:	18 (based on 220-lb. load)
Maximum angle of recline:	40°